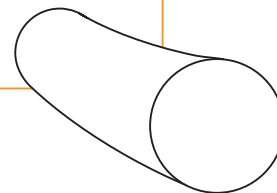




MOPYLEN® CV

- guter Knotensitz
- dauerhaft hohe Reißfestigkeit
- hervorragende Fadengleitfähigkeit
- alterungsbeständig
- SURE POINT Nadeln zur optimierten Anwendung innerhalb der cardio-vaskulären Chirurgie

- *excellent knot security*
- *consistently high tear resistance*
- *excellent tissue passage*
- *non-ageing*
- *SURE POINT needles for optimised use within cardio-vascular surgery*



Eigenschaften

- guter Knotensitz
- dauerhaft hohe Reißfestigkeit
- hervorragende Fadengleitfähigkeit
- alterungsbeständig
- SURE POINT Nadeln zur optimierten Anwendung innerhalb der cardio-vaskulären Chirurgie

MOPYLEN® CV ist ein spezieller Faden für die cardio-vaskuläre Chirurgie und wird ausschließlich mit SURE POINT Nadeln armiert.

MOPYLEN® CV ist ein synthetischer Faden, der durch Polymerisation von Propylen hergestellt wird. Die Herstellung des Fadens erfolgt aus dem gefärbten Granulat nach dem Trockenspinverfahren. Der Faden ist hydrophob, hat also praktisch keine Wasserabsorption und ist chemisch inert.

Das Material ist mit einem physiologisch unbedenklichen Farbstoff eingefärbt.

- ➔ Farbe: blau
- ➔ Chemische Bezeichnung: isotaktisches Polypropylen
- ➔ Fadenstärke: USP 8-0 – 0 (0.4-3.5 metric)
- ➔ Konfektionsangebote: Nadel-Faden-Kombinationen
- ➔ Sterilisation: Ethylenoxid

Properties

- excellent knot security
- consistently high tear resistance
- excellent tissue passage
- non-ageing
- SURE POINT needles for optimised use within cardio-vascular surgery

MOPYLEN® CV is a special suture for cardio-vascular surgery and is crimped exclusively with SURE POINT needles.

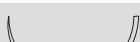
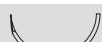
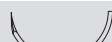
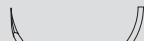
MOPYLEN® CV is a synthetic suture, which is manufactured by polymerising propylene. The suture is produced from the dyed granules using the dry spinning process. The suture is hydrophobic, i.e., it absorbs practically no water and is chemically inert.

The material is dyed with a physiologically harmless dye.

- ➔ Colour: blue
- ➔ Chemical name: isotactic polypropylene
- ➔ Thread diameter: USP 8-0 – 0 (0.4-3.5 metric)
- ➔ Types of packaging: needle-thread-combinations
- ➔ Sterilization method: ethylene oxide

MOPYLEN® CV

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	metric	Länge / length (m)	Farbe / colour	Best.-Nr. / code no.	Packungsinhalt /box content (Stück / pieces)
HR = 1/2-Kreis Rundkörper „Sure Point“ <i>HR = 1/2 circle, taper point „Sure Point“</i>						
HR 6 ● 	8-0	0.4	0,75	blau / blue	78101CVN	12
2 x HR 10 ● 	7-0	0.5	0,75	blau / blue	78102CVN	12
	6-0	0.7	0,75	blau / blue	78103CVN	12
	5-0	1	0,75	blau / blue	78104CVN	12
2 x HR 12 ● 	6-0	0.7	0,75	blau / blue	78105CVN	12
	5-0	1	0,75	blau / blue	78106CVN	12
	4-0	1.5	0,75	blau / blue	78107CVN	12
HR 17 ● 	5-0	1	0,75	blau / blue	78108CVN	12
	4-0	1.5	0,75	blau / blue	78109CVN	12
	3-0	2	0,75	blau / blue	78110CVN	12
2 x HR 17 ● 	5-0	1	0,90	blau / blue	78111CVN	12
	4-0	1.5	0,90	blau / blue	78112CVN	12
	3-0	2	0,90	blau / blue	78113CVN	12
HR 20 ● 	3-0	2	0,75	blau / blue	78114CVN	12
	2-0	3	0,75	blau / blue	78115CVN	12
HR 27 ● 	3-0	2	0,75	blau / blue	78116CVN	12
	2-0	3	0,75	blau / blue	78117CVN	12
2 x HR 27 ● 	3-0	2	0,90	blau / blue	78118CVN	12
	2-0	3	0,90	blau / blue	78119CVN	12
HRT = 1/2-Kreis schneidende Rundkörper „Sure Point“ <i>HRT = 1/2 circle, taper cutting „Sure Point“</i>						
HR 18 ▼ 	5-0	1	0,75	blau / blue	78301CVN	12
	4-0	1.5	0,75	blau / blue	78302CVN	12
	3-0	2	0,75	blau / blue	78303CVN	12
2 x HRT 18 ▼ 	5-0	1	0,90	blau / blue	78304CVN	12
	4-0	1.5	0,90	blau / blue	78305CVN	12
	3-0	2	0,90	blau / blue	78306CVN	12
HRT 20 ▼ 	3-0	2	0,75	blau / blue	78307CVN	12
	2-0	3	0,75	blau / blue	78308CVN	12
2 x HRT 20 ▼ 	3-0	2	0,90	blau / blue	78309CVN	12
	2-0	3	0,90	blau / blue	78310CVN	12
HRT 25 ▼ 	3-0	2	0,75	blau / blue	78311CVN	12
	2-0	3	0,75	blau / blue	78312CVN	12
	0	3.5	0,75	blau / blue	78313CVN	12
2 x HRT 25 ▼ 	3-0	2	0,90	blau / blue	78314CVN	12
	2-0	3	0,90	blau / blue	78315CVN	12

MOPYLEN® CV

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	metric	Länge / length (m)	Farbe / colour	Best.-Nr. / code no.	Packungsinhalt /box content (Stück / pieces)
DR = 3/8-Kreis Rundkörper „Sure Point“ <i>DR = 3/8 circle, taper point „Sure Point“</i>						
2 x DR 5 ● 	8-0	0,4	0,75	blau / blue	78001CVN	12
DR 8 ● 	7-0	0,5	0,75	blau / blue	78002CVN	12
	6-0	0,7	0,75	blau / blue	78003CVN	12
2 x DR 8 ● 	7-0	0,5	0,75	blau / blue	78004CVN	12
	6-0	0,7	0,75	blau / blue	78005CVN	12
DR 10 ● 	7-0	0,5	0,75	blau / blue	78006CVN	12
	6-0	0,7	0,75	blau / blue	78007CVN	12
	5-0	1	0,75	blau / blue	78008CVN	12
2 x DR 10 ● 	7-0	0,5	0,75	blau / blue	78009CVN	12
	6-0	0,7	0,75	blau / blue	78010CVN	12
	5-0	1	0,75	blau / blue	78011CVN	12
	6-0	0,7	0,90	blau / blue	78012CVN	12
DR 12 ● 	6-0	0,7	0,75	blau / blue	78013CVN	12
	5-0	1	0,75	blau / blue	78014CVN	12
	4-0	1,5	0,75	blau / blue	78015CVN	12
2 x DR 12 ● 	6-0	0,7	0,75	blau / blue	78016CVN	12
	5-0	1	0,75	blau / blue	78017CVN	12
	4-0	1,5	0,75	blau / blue	78018CVN	12
DR 18 ● 	5-0	1	0,75	blau / blue	78019CVN	12
	4-0	1,5	0,75	blau / blue	78020CVN	12
	3-0	2	0,75	blau / blue	78021CVN	12
2 x DR 18 ● 	5-0	1	0,90	blau / blue	78022CVN	12
	4-0	1,5	0,90	blau / blue	78023CVN	12
	3-0	2	0,90	blau / blue	78024CVN	12
DR 20 ● 	2-0	3	0,75	blau / blue	78025CVN	12
	0	3,5	0,75	blau / blue	78026CVN	12
2 x DR 20 ● 	3-0	2	0,90	blau / blue	78027CVN	12
	2-0	3	0,90	blau / blue	78028CVN	12
2 x DR 27 ● 	2-0	3	0,90	blau / blue	78029CVN	12
	0	3,5	0,90	blau / blue	78030CVN	12

MOPYLEN® CV

Sterile Nadel-Faden-Kombination / *needle suture combination, sterile*

Nadel / needle	Stärke / gauge USP	metric	Länge / length (m)	Farbe / colour	Best.-Nr. / code no.	Packungsinhalt /box content (Stück / pieces)
DRT = 3/8-Kreis schneidende Rundkörper „Sure Point“ <i>DRT = 3/8 circle, taper cutting „Sure Point“</i>						
2 x DRT 8 ▼ 	7-0	0.5	0,75	blau / blue	78201CVN	12
	6-0	0.7	0,75	blau / blue	78202CVN	12
DRT 10 ▼ 	7-0	0.5	0,75	blau / blue	78203CVN	12
	6-0	0.7	0,75	blau / blue	78204CVN	12
2 x DRT 10 ▼ 	7-0	0.5	0,75	blau / blue	78205CVN	12
	6-0	0.7	0,75	blau / blue	78206CVN	12
DRT 12 ▼ 	6-0	0.7	0,75	blau / blue	78207CVN	12
	5-0	1	0,75	blau / blue	78208CVN	12
2 x DRT 12 ▼ 	6-0	0.7	0,75	blau / blue	78209CVN	12
	5-0	1	0,75	blau / blue	78210CVN	12
	4-0	1.5	0,75	blau / blue	78211CVN	12
2 x DRT 18 ▼ 	5-0	1	0,75	blau / blue	78212CVN	12
	4-0	1.5	0,75	blau / blue	78213CVN	12
	3-0	2	0,75	blau / blue	78214CVN	12
DRT 20 ▼ 	3-0	2	0,75	blau / blue	78215CVN	12
	2-0	3	0,75	blau / blue	78216CVN	12
2 x DRT 25 ▼ 	3-0	2	0,90	blau / blue	78218CVN	12
	2-0	3	0,90	blau / blue	78219CVN	12
	0	3.5	0,90	blau / blue	78217CVN	12